


## GLP documentation

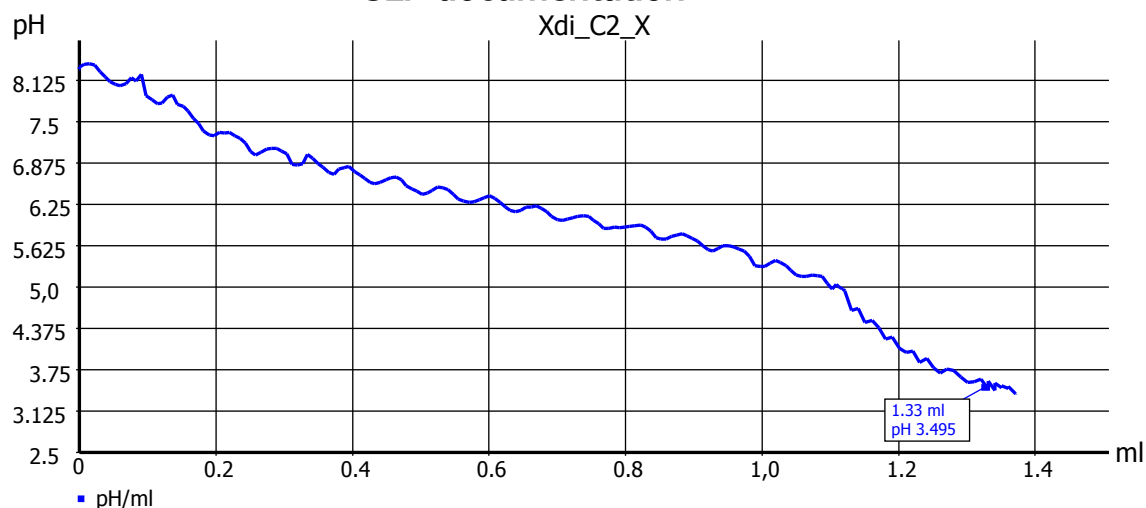

### Method data

Method name:	Alkalinity_Gran_0.1N	Titration duration:	8 m 4 s
End date:	01/25/22	End time:	14:18:36

### Titration data

Sample ID:	Xdi_C2_X	Pattern:	50.000 ml
Start pH:	pH 8.292	End pH:	pH 3.382
Start temperature:	22.2 °C (a)	End temperature:	22.4 °C (a)
Zero point:	pH 7.10 / 4.8 mV	Slope:	92.0 % / -54.4 mV/pH (a)
EP:	1.3280 ml / pH 3.495	TA:	2.677 mmol/l

### Calculation formula

TA:	$(EP1) \cdot T \cdot M \cdot F1 / (V/F2)$	Decimal places:	3
Mol (M)	1.00780		
Titre (T):	0.1 (m)	Factor 1 (F1):	1000.0000
Pattern (V):	50.0000 ml (f)	Factor 2 (F2):	1.0000
Statistics:	Off		

### Device information

Device: TitroLine 7000  
Serial number: 10012653  
Software version: 1.18.1121.247

#### Method data overall view

Method name:	Alkalinity_Gran_0.1N	Created at:	05/19/19 13:13:30
Method type:	Automatic titration	Last modification:	12/10/19 22:31:24
Measured value:	pH	Damping settings:	None
Titration mode:	End pt.	Documentation:	GLP
Linear steps:	0.0100 ml		
Measuring speed / drift:	10 s		
Initial waiting time:	0 s	Stirring titration:	free
Stirring start:	free	Titration direction:	Decrease
Pre-titration:	Off		
Pre-dosing 1:	Off	Pre-dosing 2:	Off
Endpoint	pH 3.495	delta endpoint:	pH 1.450
		Endpoint delay:	10 s

#### Dosing parameter (Titration data)

Dosing speed	2.50 % (1.00 ml/min)	Filling speed	30 s
Maximum dosing volume	50.0000 ml		

#### titration burette

Serial number:	10012653	Software version:	1.18.1121.247
Name:	TitroLine 7000		

#### Unit values

Unit size:	20 ml	Unit ID:	10043063
Reagent:	HCl 0.1M	Batch ID:	no entry
Concentration:	1.87500	Conc. determined at:	01/18/22 14:37:30
Expire date:	--	Opened/compounded:	--
Inspection according ISO:	--	Last modification:	01/18/22 14:37:35

#### Electrode data

Name:	SCP-pHT-A170MF-3M-DIN-N	Batch ID:	B210806001
Calibration:	01/25/22 9:54:44		

#### Device information

Device: TitroLine 7000  
Serial number: 10012653  
Software version: 1.18.1121.247